

ARTIA™
RECONSTRUCTIVE
TISSUE MATRIX *specifically*
designed for use in breast
plastic and reconstructive
surgeries^{1,2}

Allergan's porcine-derived
acellular dermal matrix

*Available in rectangular and contour sizes
to provide you with more options¹*

WHAT IS ARTIA™

Reconstructive Tissue Matrix?

- A **sterilized, surgical mesh that is derived from porcine skin**, and is processed and preserved in a patented aqueous phosphate-buffered solution containing matrix stabilizers²
- A **tissue patch designed to reinforce soft tissue** where weakness exists, and for the surgical repair of damaged or ruptured soft tissue membranes²
- Provides for a **strong and biocompatible implant** that will incorporate into the patient's tissue with associated cell and microvascular ingrowth²

*Product characteristics*²

- Pre-sterilized²
- Does not require refrigeration—can be stored from -8°C–30°C²
- Ready-to-use out of the package following a 2-minute soak in sterile saline or lactated Ringer's solution²
- Once opened and retained in sterile solution, must be used within 4 hours²

*Perforated options designed to:*¹

- Allow fluid to flow through matrix¹
- Support tissue ingrowth^{1,3*}
- Incorporate a 1.2 cm space around the perimeter of the perforation pattern in each piece to avoid suturing interference¹

*Contoured shape options specifically designed for breast reconstruction, and designed to:*¹

- Reduce time for product trimming¹
- Make intraoperative placement easier and more predictable¹

Consider Artia™ to **reconstruct, recontour and reform** human soft connective tissue²

* As demonstrated in a non-human primate abdominal wall model

HOW DOES ARTIA™ COMPARE to other acellular dermal matrices (ADMs)?

Artia™ enables the surgical repair of damaged or ruptured soft tissue membranes in breast plastic and reconstructive applications²


Material performance



Drape

(ratio of draped area over undraped area)

- Artia™ demonstrated **superior drapability** than Strattice™ Pliable⁴
- Artia™ demonstrated statistically similar drapability to AlloDerm™⁴

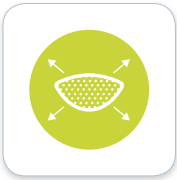

Variability in tissue stretch

(standard variability in strain)


- Artia™ demonstrated **3x less variation in stretch** than AlloDerm™⁵
- Artia™ demonstrated a **higher average strain than Strattice™ Pliable**⁵
- In terms of stretch, Artia™ behaved in a **more predictable** way than AlloDerm™ RTU⁶



Thickness uniformity


(standard deviation in thickness)

- Artia™ demonstrated a **more uniform thickness from piece to piece** than AlloDerm™⁷
- Compared to AlloDerm™, Artia™ demonstrated similar **thickness consistency within pieces**⁷

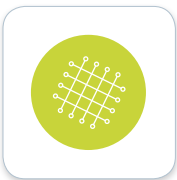


Out-of-the-package histology⁸

[hematoxylin and eosin (H&E) staining]

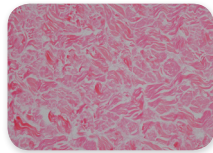

H&E staining of Artia™
Image courtesy of Allergan

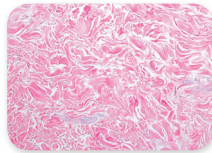

H&E staining of AlloDerm™
Image courtesy of Allergan

- Artia™ demonstrated **similar histological microstructure to AlloDerm™**, both of which are structurally **similar to human tissue**⁸

Artia™ offered equal or better handling, drapability, and ease of suturing through compared to Strattice™ Pliable and SurgiMend®^{9*}

* Based on qualitative evaluation by European surgeons (n=10).

HOW DOES ARTIA™ COMPARE to other acellular dermal matrices (ADMs)?

Preclinical results

Artia™ demonstrated superior biologic response to Strattice™ Pliable in the immune-competent-rat-subcutaneous-implant model at both 2 and 4 weeks¹⁰




Inflammatory response (Monocyte activation)


- **Human monocytes did not activate** in presence of Artia™ in vitro, similar to AlloDerm™¹¹
- Artia™ demonstrated a **minimized immunological response** in vitro¹¹



Subcutaneous in vivo testing (Immune-competent rodent, H&E staining)

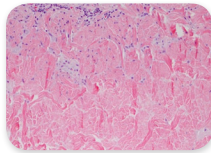

Strattice™ at 4 weeks
(10x magnification)¹²
Image courtesy of Allergan

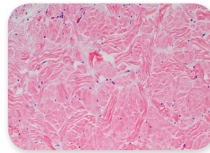

Artia™ at 2 weeks
(10x magnification)¹²
Image courtesy of Allergan

- Artia™ demonstrated **less inflammatory cell infiltration and greater and more consistent cell repopulation** (fibroblasts) compared to Strattice™¹²



Subcutaneous in vivo testing (non-human primate model)

- The overall **biologic response** of Artia™ was found to be **superior to Strattice™ Pliable**¹³
- The overall **biologic response** of Artia™ was found to be **equivalent to AlloDerm™**¹³

Artia™ is the latest addition to the ADM family and offers better biologic response compared to Strattice™ Pliable, as demonstrated in a primate model¹³



HOW DOES ARTIA™ COMPARE to other acellular dermal matrices (ADMs)?

Clinical results

In a short-term safety profile evaluating the outcomes following the use of Artia™ in consecutive immediate implant-based breast reconstruction in 17 patients,¹⁴

- There were **no implant losses** (mean follow up of 177 days), compared with over 9% identified in an audit for the Implant Breast Reconstruction evaluation (IBRA) study¹⁴
- Artia™ appeared to provide an **effective, functional, reconstructive matrix with a generally good safety profile** when utilised for either therapeutic or prophylactic mastectomy with implant reconstruction¹⁴
- Investigators stated they associated Artia™ with **minimal post-operative complications** and **good patient satisfaction**¹⁴

Outcome data in patients who had received Artia™-assisted breast reconstruction was compared to patients recipient of other ADMs such as Strattice™ or Surgimend™¹⁵

- Artia™ was **used in 31 patients** undergoing 51 implant-based breast reconstructions between July 2016 and August 2017¹⁵
- The **overall complication rate was 9.8%** with reports of complications in 5 breasts from 4 patients (mean follow-up of 171 days)¹⁵
 - Minor complications: 3 cases of seromas which were drained in the post-operative clinic appointment¹⁵
 - Major complications: 1 patient suffered bilateral implant losses following infection and skin necrosis¹⁵
- **No delays to adjuvant treatment** in therapeutic cases¹⁵
- **Complication rates compared with those of well-established ADMs** such as Strattice™, Surgimend™ and AlloDerm™ from literature¹⁵

Artia™ has been demonstrated in breast reconstruction patients to be effective, with a generally good safety profile, and with minimal post-operative complications.^{14,15}

ORDERING INFORMATION

Artia™ RTM

Allergan's porcine-derived acellular dermal matrix for breast reconstruction

Artia™ RTM'



Product Code		Product Size	Coverage
06161097EU	06161097	6 x 16 cm	96 cm ²
08161097EU	08161097	8 x 16 cm	128 cm ²
10161097EU	10161097	10 x 16 cm	160 cm ²
08201097EU	08201097	8 x 20 cm	160 cm ²

ORDERING INFORMATION

Artia™ Contour

Specifically designed for breast reconstruction¹

Artia™ Contour¹



Product Code		Product Size	Coverage
CS1097EU	CS1097	Contour Small 14,7 x 7,3	77 cm ²
CM1097EU	CM1097	Contour Medium 9.6 x 19.3 cm	132 cm ²
CL1097EU	CL1097	Contour Large 10.7 x 21.5 cm	164 cm ²
	CXL1097	Contour X-Large 23,7 x 11,8	200 cm ²

ADDITIONAL INFORMATION

Ordering:

To place an order, please call **GD Medical at 040 - 3031090**, contact your account manager or visit gdmedical.nl.

Please refer to the package insert supplied with each device for indications, contraindications, warnings, instructions for use, limited warranty, and patient, physician, and manufacturer information.

Adverse events and reporting:

Potential adverse events are those typically associated with the surgical mesh materials and/or their implantation procedures including, but not limited to, infection, foreign body response, failure to integrate, hematoma, seroma formation, lack of tissue perfusion, wound dehiscence, recurrence of tissue defect, and inflammation. If an unanticipated event occurs, alteration of surgical plan may be necessary at the surgeon's discretion.²

References:

1. Allergan Inc. Unpublished Data: LifeCell Portfolio—Product Options, INT/0218/2018. 2. Allergan Inc. Artia™ Reconstructive Tissue Matrix Instructions for Use, 2018. 3. Allergan Inc. Data on File: Evaluation of Immune/Inflammatory Response, Collagen Deposition, and Myofibroblast Presence in Perforated Porcine Acellular Dermal Matrix (pADM), INT/0260/2018. 4. Allergan Inc. Data on File: Artia™ In-Vivo testing – Draping, INT/0768/2017(1). 5. Allergan Inc. Data on File: Artia™ – Mechanical testing, INT/0770/2017. 6. Allergan Inc. Data on File: Artia™ tissue stretch, INT/0204/2018. 7. Allergan Inc. Data on File: Comparative Testing of ADMs - Thickness, INT/0225/2018. 8. Allergan Inc. Data on File: Artia™ In-Vitro Testing, out of package histology, INT/0779/2017. 9. Allergan Inc. Data on File: Artia™ Hands-on assessment, INT/0261/2018. 10. Allergan Inc. Data on File: Artia™ & Strattice™, 4 week histology, INT/0767/2017. 11. Allergan Inc. Data on File: Artia™, measure monocyte activation by flow cytometry, INT/0764/2017a. 12. Allergan Inc. Data on File: Artia™ vs. Strattice™ in an immune competent rodent - LRD2014-05-003, INT/1022/2017. 13. Allergan Inc. Data on File: Artia™ Results from non-human primate subcutaneous model, INT/0769/2017. 14. Wood B, Tait C, Linforth R. The short-term safety profile of Artia™ ADM implant based immediate breast reconstruction. Poster session presented at: Oncoplastic and Reconstructive Breast Surgery International Scientific Meeting; 2017 Sept 25-27; Nottingham UK. 2017. 15. Fakim B, Highton L, Murphy L. Breast reconstruction with Artia™. Poster session presented at: Oncoplastic and Reconstructive Breast Surgery International Scientific Meeting; 2017 Sept 25-27; Nottingham UK. 2017.

